


Bridge Breast Milk Management allows for an accurate snapshot of the status of maternal and donor milk for the newborn. Within Bridge, a summary of all active bottle activity for all patients in a location can be found. Active bottles are bottles that are not in a final status. The system displays the patient's name, room and bed, and the bottle counts in each status. The number of fresh, frozen, and thawed bottles can be viewed for each patient.

Viewing Available Patient Bottles


STEP 1: Select Milk Mgmt.



STEP 2: Select Available Patient Bottles.

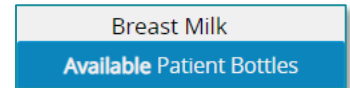

STEP 3: Use the Point of Care or Patient Status to filter the list.

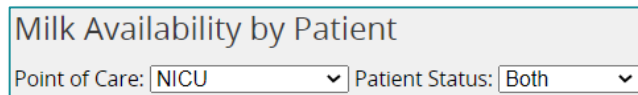

STEP 4: Select a patient and scan the patient's ID or storage label to view information.

STEP 5: If you want to dispose all bottles for the patient, select **Dispose All**. You can also dispose bottles individually by checking the box in the **Dispose** column.

Browse Bottles By Patient							
Storage State: All				Bottle #: <input type="text"/>			
Total Bottles: 4		Fresh :4 Frozen :0 Thawed :0					
Edit	Bottle # ▲	Batch # ▼	Expires ▼	State ▼	Received ▼	Print	<input type="checkbox"/> Dispose All Dispose
	1054		6/17/2021 13:30 EDT	FRESH	6/15/2021 13:30 EDT (by NICURN1) Collected Date: 6/15/2021 13:30 EDT	Print	<input type="checkbox"/>
	1055		6/17/2021 13:30 EDT	FRESH	6/15/2021 13:30 EDT (by NICURN1) Collected Date: 6/15/2021 13:30 EDT	Print	<input type="checkbox"/>

STEP 6: If donor breast milk is enabled, the **Donor Breast Milk Bottles** link is displayed. The link is available for all POCs. Select the link to view active donor bottles. If you have the privilege to correct breast milk bottles, the **Dispose** column also is available, and you can dispose bottles individually.

STEP 7: Select **Back** to dispose the selected bottles. If all bottles are disposed for a patient, the patient is no longer displayed on the list because they no longer have any remaining active bottles.