



The PCA dynamic group is located within the Lines Devices Procedures iView band. This adjustment will allow Inpatient Nursing staff to relay accurate documentation on how much medication was received via PCA for Palliative Care and Pain Management Providers to determine the appropriate oral or long-acting dosing.

PCA iView Documentation

If the PCA section is not in view, the section will need to be pulled in using **Customize View** within the Lines Devices Procedures iView band.

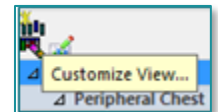


STEP 1: Create the PCA Dynamic Group by clicking the dynamic group icon within the PCA Section.

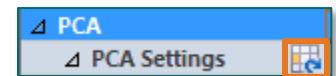

STEP 2: Select the appropriate medication from the list and click OK.

STEP 3: Enter PCA Activity by double-clicking the cell.


- Initiate and Pump Settings Change requires an RN Witness.

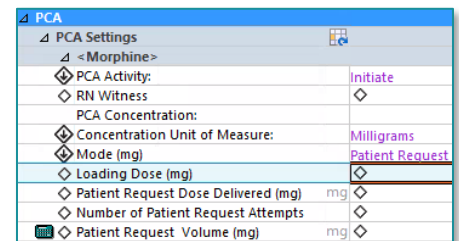

STEP 4: Enter the PCA Concentration, Concentration Unit of Measure and Mode.

NOTE: Based on the Mode selection, documentation and a calculation for Background volume, Patient Request Volume, or Background/Patient Request Volume populates. If any documentation field under Mode are left blank, the formula will not generate the calculation. For example: *If the patient does not receive a Loading dose, a zero (0) must be entered.* Hover over the calculation icon to view the formulas.

STEP 5: Document the PCA fluid volume in the IV Intake section of the Intake and Output band. If not in view, click **Customize** and select **Other Medication 1** or **Other Medication 2**, then OK.

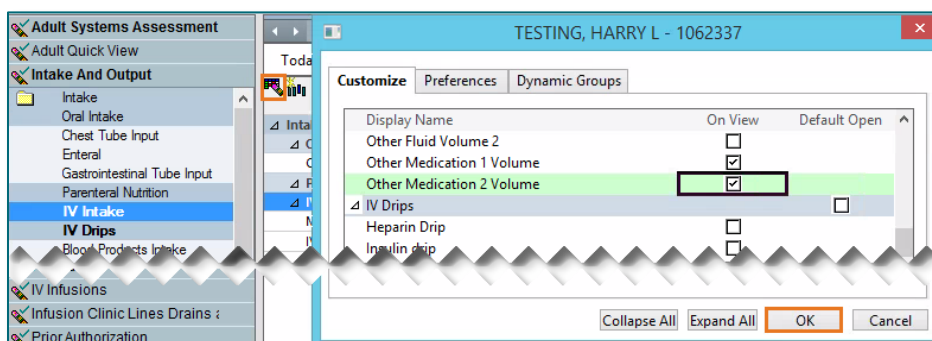

➤ The Alertness Score for the PCA checks will be documented in the **Pain Assessment** section of the applicable **System Assessment** iView band.

